




Tip Top™ Dispenser Caps

A safe and cost effective way of dispensing samples in your lab.


The Tip Top™ is a disposable cone shaped cap that fits all standard 13mm or 16mm blood collection tubes. A small opening allows you to dispense a sample without removing the Tip Top™ cap so you can prepare smears on slides or dispense into reagent wells or microcentrifuge tubes safely.

- **Convenient** - just tilt the tube and squeeze to dispense one drop
- **Safe** - you don't need to open the tube and pipet out of the tube so spillage is reduced along with potential exposure to infectious disease
- **Low Evaporation** - A flexible cap is available to cover the Tip Top™ so you can store samples in the refrigerator.

The Tip Top™ cap can be used with almost any liquid. Made from chemically resistant polyethylene they can dispense serum, reagents, plasma, whole blood and Cerebrospinal fluid. The integral porous polyethylene filter efficiently separates serum and can also be used as a reagent blender.

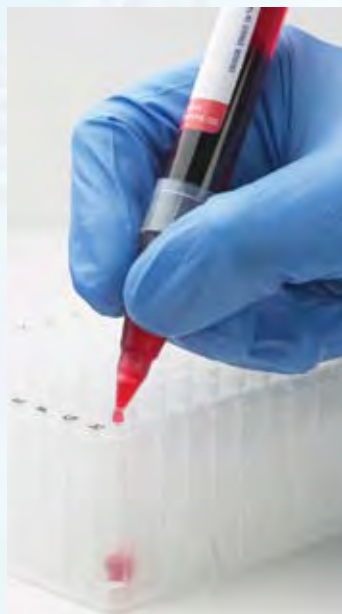

- Fits standard 13 x 100 and 16 x 100 collection tubes
- With integral filters for prep work

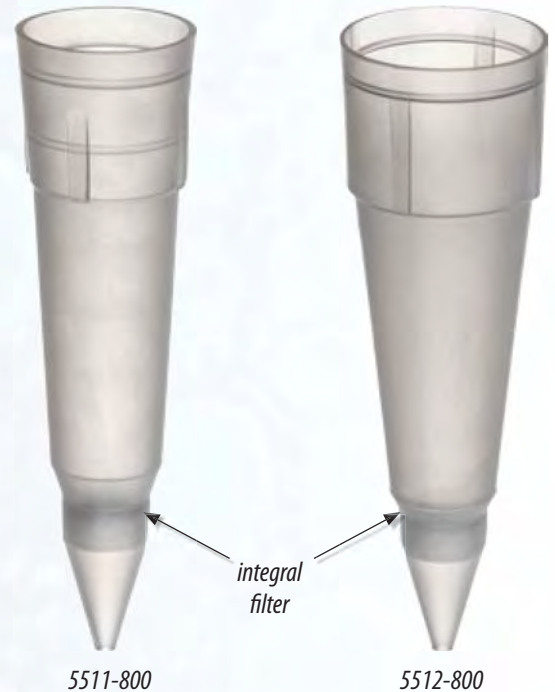

Tip Top™ Dispenser Caps

Item	Description	Qty/Pk	Qty/Cs
5511-800	Natural, Fits 13mm tubes, 1000/bag	1000	5000
5512-800	Natural, Fits 16mm tubes, 1000/bag	1000	5000

Closures for Sealing Tip Top™ Caps

Item	Description	Qty/Pk	Qty/Cs
5515-804	Closure, Red, 1000/bag	1000	5000